

Features built for mobility



VeinViewer® Flex packs the latest advances in superior vascular imaging into a handheld, mobile package. VeinViewer can help deliver faster IV procedures on your difficult patients through greatly improved first stick success.



Fast Swap Flexible Power

Flex runs on Fast Swap Battery and AC power so there are no concerns about a down device. Batteries are compatible with Vision2 for system-wide optimization.

S-Mount (Stability Mount)

When you need to be more hands-on, the attachable S-mount gives you a hands-free approach for an Eyes On Patient™ technique.

Hospital Grade Durability

Durability withstands even the busiest shifts in the most demanding conditions.

Flex Imaging Unit

Built to be controlled and positioned easily with one hand, Flex is ergonomically designed for simplified procedures. The standard Universal Mode provides HD and Df² projection technology, a clear and accurate view of superficial vein patterns in real time, even in remote locations. Gain even more customizable features with the ASSESS™ Imaging Suite, available with a licensed upgrade.



Standard and optional accessories shown.

VeinViewer Flex Available Accessories


- S 2 S-Mounts
 - Super Clamp
 - C-Clamp
- S 2 Fast Swap Batteries
- S Battery Charger
- S AC Power Adapter
- S International Power Adapters
- S Nylon Storage Case
- S USB Cable
- S Product Training DVD
- O Travel Case
- O Car Charger

S - Standard O - Optional




Every day you have more and more patients with challenging IV start needs. VeinViewer provides a real solution to help meet those needs.


Flex gives you revolutionary HD views in real-time, so you have the best view of vasculature procedures Pre-, During- and Post-access. MaxDepth allows a view of peripheral veins up to 10mm deep and blood patterns up to 15mm beneath the surface of the skin.

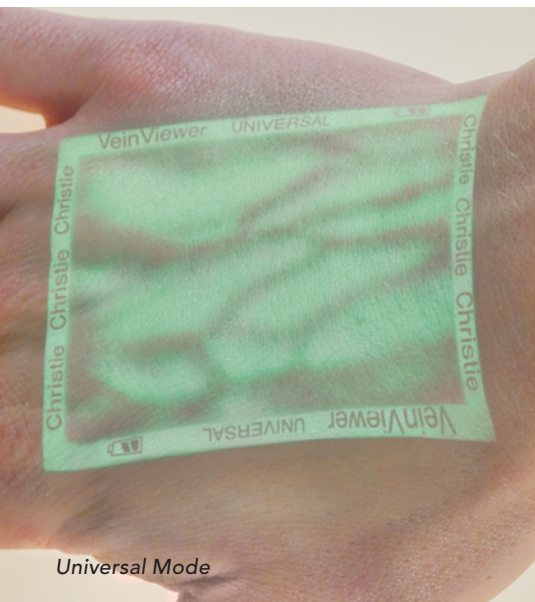

Universal Mode

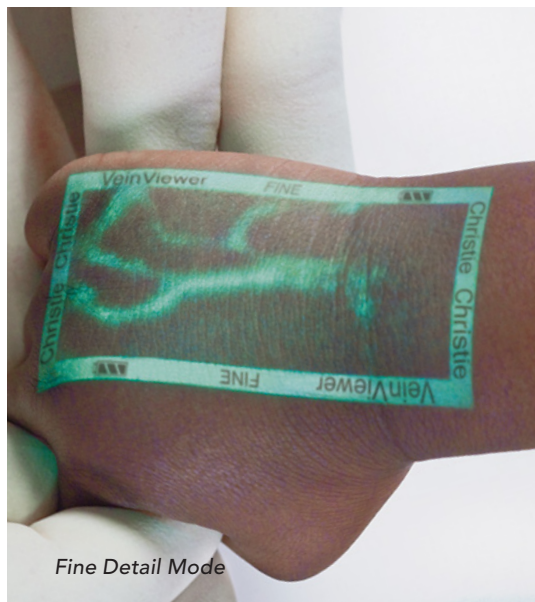

Fine Detail Mode

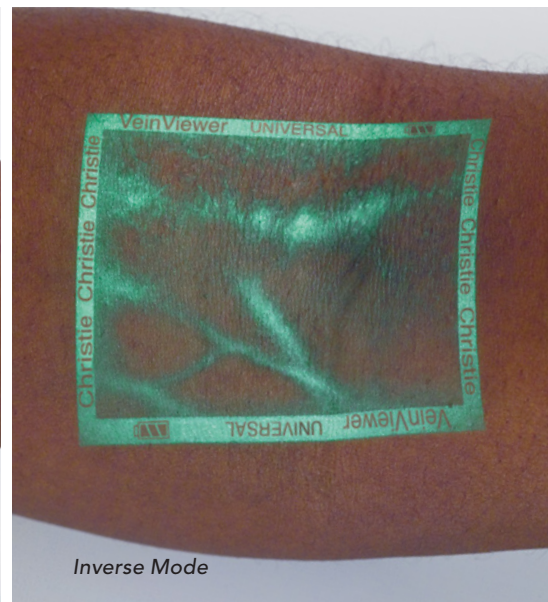

Inverse Mode

ASSESS™ Imaging Suite

See veins like never before with the ASSESS Imaging Suite of software only available on VeinViewer products. The customizable modes let you choose the best view for each patient; giving you more control to ensure a better assessment for every PIV procedure.

ASSESS™

- **Universal** - useful for all patients
- **Fine Detail*** - critical for seeing small veins in pediatric or sclerotherapy patients
- **Inverse*** - reverses the color of the veins within the image field for further image customization
- **Resize*** - choose 1 of 3 window sizes, a great feature for pediatric patients
- **MaxBright**** - brightest vein illumination image, needed for varying lighting conditions
- **Image Capture*** - HIPAA compliant, takes up to 200 static .PNG images of the blood pattern for export to EMR documentation

*Requires software upgrade.

** Requires software upgrade. Available in the EMS configuration only.

Your skills. Our technology. Better patient outcomes.

info@veinviewer.com

www.veinviewer.com

1-877-SEE-VEIN (1-877-733-8346) TOLL-FREE